TOOTH DESENSITIZER

TEETHMATE™ DESENSITIZER

CE0197

ENGLISH INSTRUCTIONS FOR USE

- I. INTRODUCTION
- TEETHMATE DESENSITIZER is composed of POWDER which contains Calcium Phosphate as a major component and LIQUID which mainly contains water. POWDER and LIQUID are mixed to obtain slurry. The slurry applied to the affected area transforms to a hardened material of hydroxyapatite and it suppresses hypersensitivity by sealing the dentin tubules and microcracks in the enamel.
- II. INDICATIONS
- TEETHMATE DESENSITIZER is indicated for reduction of tooth hypersensitivity by the following treatments:
- [1] Treatment of dentin exposed by toothbrush abrasion, gingival recession, periodontal disease and/or acid erosion

[2] Treatment of dentin after mechanical tooth cleaning, scaling and/or root planing

[3] Treatment of tooth surface before and/or after bleaching

[4] Treatment of prepared dentin for fillings and/or prosthetic restorations
- III. CONTRAINDICATIONS
- Patients with a history of hypersensitivity to this product
- IV. PRECAUTIONS
1. Safety precautions

1. This product contains substances that may cause allergic reactions. Avoid use of the product in patients with known allergies to any of the components.

2. If the patient demonstrates a hypersensitivity reaction, such as rash, eczema, features of inflammation, ulcer, swelling, itching or numbness, discontinue use of the product and seek medical attention.

3. Avoid direct contact with the skin and/or soft tissue to prevent hypersensitivity. Wear gloves or take appropriate precautions when using the product.

4. Use caution to prevent the product from coming in contact with the skin or getting into the eye. Before using the product, cover the patient's eyes with a towel to protect them in the event of splashing material.

5. If the product comes in contact with human body tissues, take the following actions:
<If the product gets in the eye>
Immediately wash the eye with copious amounts of water and consult a physician.
<If the product comes in contact with the skin or the oral mucosa>
Immediately wipe the area with a cotton pellet or a gauze pad moistened with alcohol, and rinse with copious amounts of water.

6. Use caution to prevent the patient from accidentally swallowing the product.

7. Do not reuse the applicator brush and the interdental brush to avoid cross contamination. Discard them after use.

8. The patient should be instructed to avoid drinking or eating for 45 minutes after the completion of treatment with the product. Or the sealing material constructed by the product may drop out from the dentin tubules and the sealing effect may decrease.

9. This product does not permanently suppress hypersensitivity. If the hypersensitivity recurs, repeat treatment with this product.

2. Handling and manipulation precautions

1. The product must not be used for any purposes other than specified in [II.INDICATIONS].

2. The use of this product is restricted to licensed dental professionals.

3. Use a pulp capping agent in a cavity close to the pulp.

4. Do not mix the product with any other dental materials.

5. If the instrument or the container of this product is damaged, take care of yourself from any damages and immediately discontinue use.

6. Shake the powder container gently before measuring. Measure a level spoonful of POWDER by using the inner cap edge.

7. After dispensing LIQUID or measuring POWDER, each container should be stoppered tightly as soon as possible. In particular, in the powder container, take care not to allow moisture infiltrating because the powder component reacts to water to harden.

8. If the hardened material is found in the powder container due to reacted to water, remove the hardened part before use.

9. When dispensing the liquid component into the mixing dish, hold the container completely vertically. Or the correct amount of LIQUID is not dispensed.

3. Storage precautions

1. The product must be used by the expiration date indicated on the package.

2. The product must be stored at 2 - 25°C / 36 - 77°F when not in use.
When the powder container is stored in the refrigerator, it should be left unopened at room temperature for more than 15 minutes after taking out from the refrigerator, before each use. Less exposure time at room temperature may cause condensation on the inside of the powder container.

3. The product must be kept away from extreme heat or direct sunlight.

4. The product must be stored in a proper place where only dental practitioners can access.

V. COMPONENTS

Please see the outside of the package for contents and quantity.

Principal ingredients

1) POWDER

• Tetracalcium phosphate

• Dicalcium phosphate anhydrous

2) LIQUID

• Water

• Preservative

3) Accessories

• Mixing dish

• Measuring spoon

• Applicator brush (superfine <white>)

VI. CLINICAL PROCEDURES

A. Standard procedure I (Indications [1] to [3])

[1] Treatment of dentin exposed by toothbrush abrasion, gingival recession, periodontal disease and/or acid erosion

[2] Treatment of dentin after mechanical tooth cleaning, scaling and/or root planing

[3] Treatment of tooth surface before and/or after bleaching

A-1. Tooth surface cleaning

Clean the tooth surface in the usual manner. Wipe away saliva and plaque from the affected area using cotton pellets etc. If the affected area is severely contaminated with plaque or stain, use a suitable tooth surface abrasive. After tooth surface cleaning, the drying of the affected area by air blowing is not necessary.

A-2. Mixing procedure


Standard mixing ratio

Objective area	LIQUID	POWDER
1 tooth	1 drop	1 spoonful (small measuring spoon)
2 or 3 teeth	2 drops	1 spoonful (large measuring spoon)

(1) Hold the liquid container vertically over the mixing dish and dispense the intended drops of LIQUID into the dish.

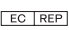(2) Shake the powder container gently before measuring. Using the measuring spoon, measure a level spoonful of POWDER by using the inner cap edge. And put it on the formerly dispensed LIQUID in the mixing dish.
- (3) Mix LIQUID and POWDER for more than 15 seconds to make the slurry, using the applicator brush.
- A-3. Application of the slurry**
Apply the slurry to the affected area using the applicator brush and then let it penetrate sufficiently by rubbing the area for more than 30 seconds. A cotton pellet can be used instead of applicator brush for large affected areas. If the applicator brush does not work properly in interdental spaces, an interdental brush can be used for this purpose.
- A-4. Rinsing**
Rinse the excess slurry with water spray or by having the patient rinse.
- A-5. Checking of the desensitizing effect**
Check the effect of treatment by blowing air onto the affected area.
* If dentin hypersensitivity persists, repeat steps A-3, A-4 and A-5.
- B. Standard procedure II (Indications [4])**
[4] Treatment of prepared dentin for fillings and/or prosthetic restorations
- B-1. Preparing cavity or abutment tooth**
Using the usual manner, prepare, clean and disinfect the cavity or abutment tooth.
- B-2. Mixing procedure**
Use the same procedure described in "A-2".
- B-3. Application of the slurry**
Use the same procedure described in "A-3".
- B-4. Rinsing**
Use the same procedure described in "A-4".
- B-5. Checking of the desensitizing effect**
Use the same procedure described in "A-5".
- B-6. Cleaning prepared dentin**
Clean the prepared dentin for more than 10 seconds using cotton pellets moistened with water to remove the minute remaining particles of the slurry from the dentin surface.
- B-7. Application of dental adhesive for filling material or cementation material for prosthetic restoration**
Follow the manufacturer's recommended procedure for each dental adhesive or cementation material.
- Note: When using dental adhesives or cementation materials which require prior phosphoric acid etching (ex: CLEARFIL PHOTO BOND), the etching step should be performed in the usual manner before applying this product.
- [WARRANTY]
Kuraray Noritake Dental Inc. will replace any product that is proven to be defective. Kuraray Noritake Dental Inc. does not accept liability for any loss or damage, direct, consequential or special, arising out of the application or use of or the inability to use these products. Before using, the user shall determine the suitability of the products for the intended use and the user assumes all risk and liability whatsoever in connection therewith.
- [NOTE]
"CLEARFIL" and "TEETHMATE" are trademarks of KURARAY CO.,LTD.
- 

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" CLEARFIL " e " TEETHMATE " sono marchi di KURARAY CO., LTD.
Kuraray Noritake Dental Inc. 1621 Sakazu, Kurashiki, Okayama 710-0801, Japan